

neoV 980 and neoV 1470 for ENT

Surgical power was never so small



The neoV surgical laser platform with your choice of 980nm or 1470nm - The smallest, lightest weight, and most portable laser system of its power class.



DESIGNED FOR SURGERY





The neoV Laser

The World's Smallest 25 Watt Laser

The neoLaser team, with over 20 years of combined medical laser industry experience, and having brought to successful launch pioneering laser technology, is now embarking on a new and

exciting path of innovation.

neoLaser is proud to launch the neoV laser platform, with its groundbreaking

design. With the ultimate

commitment to innovtive design,

detail oriented craftsmanship, and top-notch quality, neoLaser and its products will stand apart from the industry. We are here to serve your surgical needs and to work with you to develop the optimal and best designed laser based surgical solutions for your practice. We welcome your business.

Gil Shapira, Founder President & CEO



Good things come in small packages and the neoV laser is no exception.

This small, portable, superbly designed unit delivers an accurate 1470nm laser beam or a powerful 25W 980nm energy for your surgical applications in ENT.



KEY FEATURES

980 or 1470 wavelength

Up to 28W at 980 or 12W at 1470 from laser head

300 μm to 600 μm fibers

Indications in ENT

CORONA ENT/OTO fibers

EASE OF USE

Large high-brightness color touch screen

Intuitive user interface

Plug and play fiber connection

Rapid setup with surgeon presets

UNIQUE VALUE

Enhanced safety features


High reliability and low maintenance

Portability

Standard power connection



The interaction/wavelength choice

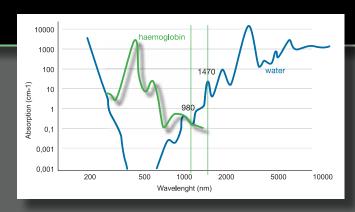

Laser interaction with soft tissue depends on the wavelength of choice.

The neoV laser comes in 2 versions – the 980nm and the 1470nm.

The 980nm laser with equal absorption in water and blood, offers a powerful all-purpose surgical tool for incision, excision, coagulation and vaporization, as well as a high power source for rapid surgical work.

The 1470nm laser with significantly higher absorption in water, provides a precision instrument for fine surgical work with reduced collateral damage around critical structures. Whatever your surgical needs, the neoV platform can deliver what you need.

Tissue interaction can be easily manipulated by changing the distance from the tissue, from precise incision when in contact mode, to coagulation and ablation when removed from the tissue.





Wide array of surgical interventions



The neoV laser system enables a wide array of surgical interventions with precision and safety. Choose either a 980nm or 1470nm system depending on the type of applications you seek. The system can be easily used in the OR or in a clinic setting, and enables a range of power levels, from delicate work near critical structures, to rapid removal of tissue.

Potential applications:

- Middle ear surgery,
 stapedotomy
- Laryngeal cancer and coredctomy
- Orapharyngeal lesions, tumor vaporization
- Tonsillectomy

- Glossectomy
- Nasal septum, polyps, cysts, mucoceles
- Subglottic stenosis and lesions
- Tracheal stenosis and lesions
- Laser assisted
 Uvulopalatoplasty






Otology - Precision and Safety



The neoV laser offers a delicate and high precision modality for inner ear treatments. By use of a special design small core fiber disposable kit, fine ablation and tissue removal can be achieved. The use of a hand-held fiber delivery system with a clearly visible aiming beam ensures safe application and high precision control of energy application.

neoLaser offers a flexible reusable hand piece system consisting of a rigid ergonomic handpiece and various termination cannulas, or a fully disposable sterile handpiece with an integrated high precision fiber.

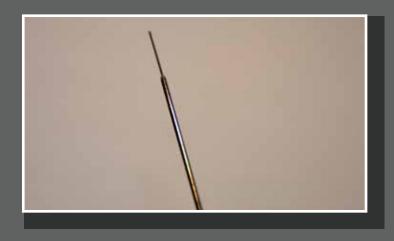





Fiber Delivery and Accessories



neoLaser offers a wide and comprehensive set of surgical accessories for multiple specialties. neoLaser's modular surgical hand piece system allows you to pick and choose the size and shape of cannula matching the treated anatomy. Simply plug your cannula of choice on the ergonomic hand piece for easy access. Custom tubes can be designed upon request.

Additional specialty hand pieces are also available, such as the laryngeal hand pieces and the suction hand piece.



Safety accessories are also available.





DESIGNED FOR SURGERY



DANGER

Av. João Cabral de Mello Neto, 850 sl:222 Barra da Tijuca - Rio de Janeiro - RJ

Tel: +55 21 2289-1210 E-mail: dp@ecmax.com.br

Proprietary

CW, Pulsed

19 VDC, 4.7 A

10 x 22 x 22 (cm)

Fiber Connection

Dimensions (HxWxD)

Operating Modes Power Requirements

VISIBLE AND INVISIBLE LASER RADIATION AVOID EYE OR SKIN EXPOSURE TO DIRECT OR SCATTERED RADIATION LASER CLASS 4 EN 60825-1-2007 P₀ = 28 Watt, 1 = 980-1470 nm P₀ < 5 mW, 1 = 532nm / 650nm